



Phonak ComPilot and BlackBerry Bluetooth pairing guide


PHONAK

life is on

1. Charge and turn on the mobile phone and ComPilot.
2. Place them next to each other.

3. Prepare ComPilot for pairing

- Press and hold the Connect  and Volume  buttons simultaneously for 2 seconds until the Audio indicator starts to rapidly blink blue.
- Start the pairing process on the phone as described below.
- Please consult your phone's user guide if you need more information.

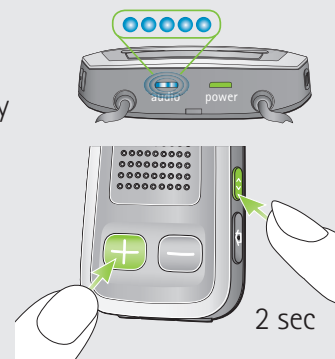

4. Prepare the mobile phone for pairing

Find the connectivity settings of the phone. With most models, the following navigation applies. The **Bold** terms might be different depending on the exact phone model.


Pairing older models

- Highlight and select **All > Options > Networks and Connections > Bluetooth Connections**
- Ensure Bluetooth is turned **On**
- Select **Add New Device > Search**
- Wait until **Phonak ComPilot** is found
- Select the **ComPilot**
- Press **Select**
- Enter the PIN code **0000** (4 zeros)
- Select **Yes**

Pairing smartphone models (touchscreen)

- Scroll to and touch **Options > Networks and Connections > Bluetooth Connections**
- Ensure Bluetooth is turned **On**
- Touch **Add New Device > Search**
- Wait until **Phonak ComPilot** is found
- Touch **ComPilot**
- Enter the passkey **0000** (4 zeros)
- Touch **Enter > Yes**

5. End of pairing

- After successful pairing, the ComPilot Audio indicator will stop rapid flashing blue.
- Additionally a headset symbol  may appear on the main screen of your phone.
- After 2 minutes the ComPilot will automatically exit from pairing mode. If step 4 was not completed before, simply start the pairing mode on ComPilot again as described in step 3.