

Audéo MINI

Product information



Small and easy, Audéo MINI is designed for hearing instrument users with mild to severe hearing losses. Audéo MINI offers industry leading CORE Performance features

including SoundRecover. Audéo MINI is available in three performance/price levels, has a choice of standard or power receiver and a range of dome and tip options.

Features/Functions

	Audéo MINI IX	Audéo MINI V	Audéo MINI III
SoundFlow	Premium	Advanced	Standard
WhistleBlock Technology	■	■	■
SoundRecover	■	■	■
Directionality	VoiceZoom (33 channels)	VoiceZoom (12 channels)	digital AudioZoom
Real Ear Sound	■	■	
SoundRelax	■		
EchoBlock System	■		
WindBlock Management	■	■	
NoiseBlock Processing	20 channels	16 channels	6 channels
DataLogging	■	■	■

Fitting software


- iPGF version 2.5 or higher, NOAH compatible

Hardware


- PC (IBM compatible) with NOAHlink and/or HI-PRO interface

Programming

Connecting the programming cable to an Audéo MINI

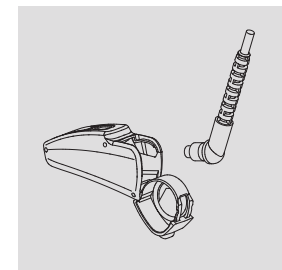

1. Open the battery compartment.

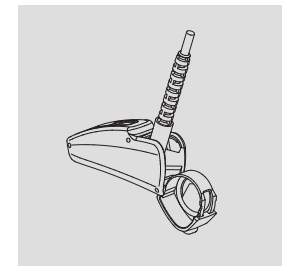

2. Plug the programming plug into the socket.
Use the programming cable CS-44 A (Phonak version) only.

