


Phonak

Roger™ Table Mic II iN

Roger Table Mic II iN is a microphone dedicated for working adults who participate in various meetings. It selects the person who's talking and switches automatically. Multiple table microphones can be connected to create a network, making it ideal for large meeting configurations. It can also transmit the sound of multi-media e.g. computer.

Technical Data

Type	Roger microphone
Dimensions (Ø x H)	69 x 69 x 14 mm (2.7 x 2.7 x 0.6")
Weight	103 g (3.63 oz)
Color options	Champagne
Included accessories	Power supply with USB charging cable, audio cable, Cinch adapter, soft case
Operating conditions	0° to +45° Celsius (+32° to +113° Fahrenheit) Relative humidity of <90% (non-condensing)
Transport and storage conditions	-20° to +60° Celsius (-4° to +140° Fahrenheit) Relative humidity of <95% for a long period of time

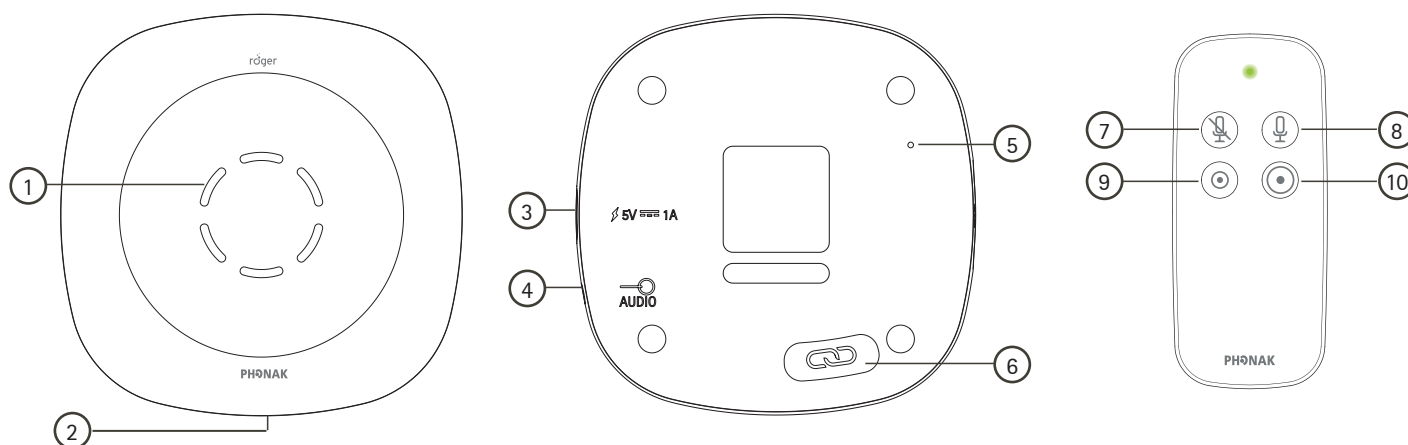

Device description

1	Microphones and indicator lights
2	On/off button
3	Micro-USB socket
4	Audio input
5	Install button (for HCP only)

6	Connect button
7	Mute on
8	Mute off
9	Short pick-up range
10	Long pick-up range



A Sonova brand

PHONAK
life is on

Features

- MultiBeam Technology for improved speech in noise
- MultiTalker Network for easy coverage of large meetings
- Remote control for muting and adjusting pick-up range
- Audio input for connection to multimedia devices
- 16 hours operating time

Roger installation

Roger Table Mic II iN contains two Roger receivers that can be installed into a pair (two) of hearing devices with RogerDirect™ without attaching an external receiver.

Microphone characteristics

Number of microphones	3
Microphone modes	Table mode

Radio information

Streaming technology	Roger technology
Frequency range	2.4 GHz ISM band including adaptive frequency hopping
RF power	70 mW max.
Transmission range (clear line of sight, free field)	Up to 80 meters / 260 feet (to hearing device with RogerDirect) Up to 40 meters / 130 feet (to Roger receiver)

Audio information

Audio input	Analog
Audio bandwidth	100 Hz – 7.3 kHz

Power supply

Type	DC, stabilized
Primary voltage	100-240 VAC
Secondary voltage	5 VDC
Max. current	1.0 A
Plug	micro-USB
Important notice	Only use the original USB Phonak power supply 5VDC, > 500 mA

Battery

Type	Lithium Polymer
Capacity	1000 mAh
Voltage	3.7 V
Operating time	16 hours

Receiver compatibility

Compatible Roger receivers	All personal Roger receivers, except SoundField
Max. connected receivers	Unlimited

Microphone compatibility (MultiTalker Network)

Compatible Roger microphones	Roger Table Mic II, Roger Table Mic II iN, Roger Select™, Roger Select™ iN, Roger Pen™, Roger Pen™ iN, Roger EasyPen and Roger Clip-On Mic
Max. connected microphones	10

Standards

Europe	EN 300 328, EN 301 489, EN 60950-1, IEC/EN 62479
Canada	RSS-247, RSS-102
Japan	ARIB-T66
USA	CFR 47, part 15.231, part 15.247, FCC OET Bulletin No.65