



Phonak Naída™ L ●

Product information

Phonak Naída Lumity is the comprehensive BTE portfolio for mild to profound hearing losses.

Powered by Phonak SmartSpeech™ Technology, Naída L hearing aids are the choice for clients who need powerful solutions on our most advanced platform Phonak Lumity.





Product description



Naída L-UP
Naída L-UP Trial
L90/L70/L50/L30



Naída L-SP
Naída L-SP Trial
L90/L70/L50/L30




Naída L-PR
Naída L-PR Trial
L90/L70/L50/L30

Battery size		675	13	built-in rechargeable
Telecoil		•	•	
Compatibility with CROS L				•
Bluetooth®		•	•	•
Hydrophobic coating		•	•	•
IP rating		IP68 ¹	IP68 ¹	IP68 ¹
Dimensions	L x W x D	42.6 x 18.8 x 8.6 mm (1.68 x 0.74 x 0.34")	40.0 x 16.9 x 8.5 mm (1.57 x 0.67 x 0.33")	36.7 x 16.4 x 8.5 mm (1.44 x 0.64 x 0.33")
	Weight	4.60 g / 0.162 oz	4.40 g / 0.155 oz	3.68 g / 0.130 oz
Charging times ²	0% to 100%			3 h
	0% to 80%			1.5 h
	30 min. fitting			10 min.

		HE11	HE11 680	HE11	HE11 680	SlimTube 4.0	HE11 680	SlimTube 4.0
Max. Power Output (dB SPL)	2 cc coupler	141	134	139	133	135	130	129
Max. gain (dB)	2 cc coupler	84	78	81	75	70	68	69
Frequency range (Hz)		< 100 - > 4700	< 100 - > 4700	< 100 - > 4900	< 100 - > 5700	< 100 - > 7100	< 100 - > 5700	< 100 - > 7100

Fitting range

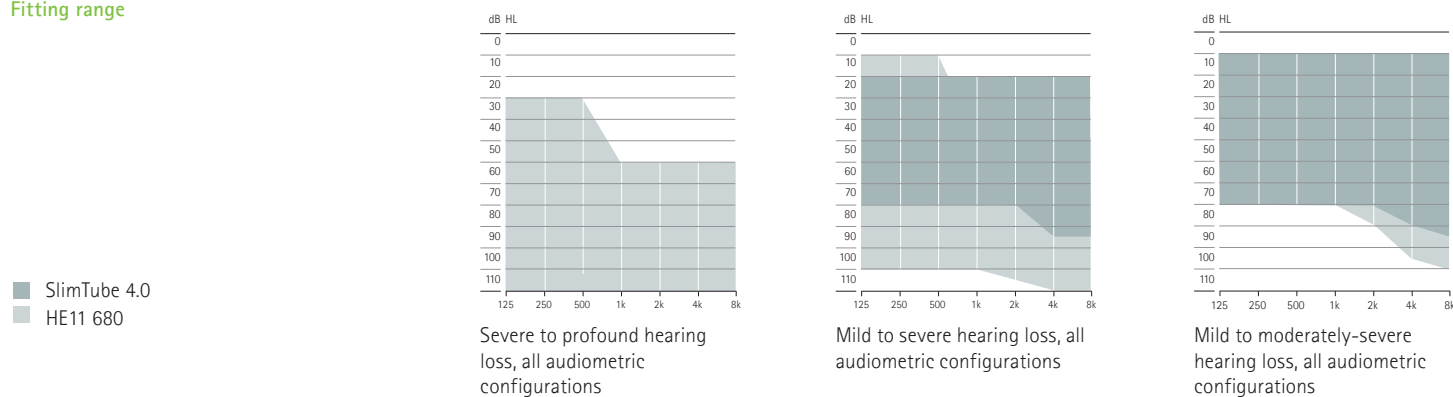

Housing colors



H0
Beige



P1
Sand Beige



P3
Sandalwood



P4
Chestnut



P5
Champagne



P6
Silver Gray



P7
Graphite Gray



P8
Velvet Black

¹ IP68 indicates that the hearing aid is water and dust resistant. It survived continuous immersion in 1 meter of fresh water for 60 minutes and 8 hours in a dust chamber as per the IEC60529 standard.
² Actual charging time depends on the remaining battery life and may be affected by temperature conditions.

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Sonova AG is under license.





Performance levels

AutoSense OS™ 5.0	Premium L90	Advanced L70	Standard L50	Essential L30
Speech in Loud Noise	•			
Speech in Car	•			
Comfort in Echo	•			
Music	•	•		
Comfort in Noise	•	•	•	
Speech in Noise	•	•	•	•
Calm Situation	•	•	•	•
Media speech + mic	•	•	•	•
Media music + mic	•	•	•	•
Additional programs				
Max. additional programs	4	4	4	4
Comfort in Echo	•			
Speech in Loud Noise	•	•		
Speech in 360°	•	•		
Speech in Noise	•	•	•	•
Calm Situation	•	•	•	•
Comfort in Noise	•	•	•	•
Music	•	•	•	•
Acoustic phone	•	•	•	•
Custom program	•	•	•	•
Phone via T-coil ± mic¹	•	•	•	•
Public T-coil ± mic¹	•	•	•	•
Streaming programs				
RogerDirect + mic	•	•	•	•
Phonak PartnerMic™ + mic	•	•	•	•
Phone call + mic	•	•	•	•

SmartSpeech Technology

Features	Premium L90	Advanced L70	Standard L50	Essential L30
SpeechSensor	•			
Speech Enhancer	•			
StereoZoom 2.0	•	•		
Dynamic Noise Cancellation	•	•		
Motion Sensor Hearing²	•	•	•	
RogerDirect™	•	•	•	•
EchoBlock	•			
Tap control²				
- Siri/Google Assistant	•	•		
- Pause/resume streaming	•	•		
- Accept/end phone call	•	•		
SoundRelax	•	•		
DuoPhone	•	•	•	
Real Ear Sound	•	•	•	
UltraZoom	•	•	•	•
Roger™ and Directional setting	•	•	•	•
Broadband Booster³	•	•	•	•
Health functionalities²⁴	•	•	•	•
WindBlock	•	•	•	•
User Preference Tuning	•	•	•	•
Tinnitus Balance⁵	•	•	•	•
Auto Acclimatization	•	•	•	•
Environmental balance	•	•	•	•
SoundRecover2	•	•	•	•
WhistleBlock	•	•	•	•
NoiseBlock	•	•	•	•
QuickSync	•	•	•	•
AOV	•	•	•	•
Fine tuning channels	20	20	16	12

1 Not available for Phonak Naída L-PR
2 Available for Phonak Naída L-PR only
3 Available for Phonak Naída L-UP only
4 Available in combination with myPhonak
5 Not available for Phonak Naída L-UP

Fitting

	Naída L-UP Naída L-UP Trial L90/L70/L50/L30	Naída L-SP Naída L-SP Trial L90/L70/L50/L30	Naída L-PR Naída L-PR Trial L90/L70/L50/L30
Phonak Target	9.0 or higher	9.1 or higher	9.0 or higher
Noahlink Wireless / Noahlink Wireless 2	•	•	•
NOAHLINK	•	•	
HI-PRO, HI-PRO2	•	•	





Acoustic coupling

			Naída L-UP Naída L-UP Trial L90/L70/L50/L30	Naída L-SP Naída L-SP Trial L90/L70/L50/L30	Naída L-PR Naída L-PR Trial L90/L70/L50/L30
Sound outlet	Hook	Damped (HE11 680), transparent	•	•	•
	Hook	Undamped (HE11), transparent	•	•	
	SlimTube 4.0	Available in left and right in five different lengths (00,0,1,2,3)		•	•
Earpiece for SlimTube	Dome	Cap dome		•	•
		Open dome Available in three sizes (S, M, L)		•	•
		Vented dome Available in three sizes (S, M, L)		•	•
		Power dome Available in three sizes (S, M, L)		•	•
	SlimTip	Acrylic		•	•
		Silicone		•	•
Earpiece for hook	Earmold		•	•	•

Phonak Naída L Trial

The Phonak Naída L Trial hearing aids can be quickly set to the required performance level in order to suit your recommendation. This provides you the ability to fit your client with the optimal hearing solution already at the first appointment or to provide a replacement when the client's hearing aid is in for repair. The Naída L Trial devices also include a pre-installed Roger receiver for easy demonstration and trial of any Roger microphone.



Phonak Naída L-PR charging options

Phonak Charger Combi BTE / Phonak Charger Combi BTE 2

Naída L-PR has all the Lumity features you love as well as the benefit of rechargeability. Clients can enjoy a full day of hearing including streaming*.

* Expected battery life of 16 hours, including 4 hours of BT classic streaming and 4 hours of TV Connector streaming.



Phonak Naída L Bluetooth® versions

Direct calls/audio streaming

Bluetooth 4.2 wireless technology and most older Bluetooth phones

Wireless communication portfolio

	Naída L-UP Naída L-UP Trial L90/L70/L50/L30	Naída L-SP Naída L-SP Trial L90/L70/L50/L30	Naída L-PR Naída L-PR Trial L90/L70/L50/L30
Wireless accessories			
TV Connector	•	•	•
Phonak PartnerMic	•	•	•
Phonak RemoteControl	•	•	•
Roger			
Roger microphones	•	•	•
Apps			
myPhonak	•	•	•
myPhonak Junior	•	•	•

For more information please visit www.phonak.com/naida-lumity

For information on Roger microphones please visit www.phonak.com/roger

For more information on myPhonak and myPhonak Junior please visit www.phonak.com/apps





life is on


Since 1947, Phonak is dedicated to preserving life quality by opening new acoustic worlds. Already back then, in the days of our foundation, our company was driven by a formative conviction: We believe that well-hearing equates to well-being and thus is essential for living life to the fullest. In fact, the sense of hearing is directly linked to social, emotional, cognitive and physical well-being. Today as in future, we thrive to offer the broadest portfolio of innovative hearing solutions. And, together with our hearing care professionals, we keep on focusing on what matters most: improving speech understanding, changing people's lives and having a positive effect on society as a whole.

www.phonak.com



027-0699-02/V3.00/2024-02/na © 2024 Sonova AG All rights reserved

sonova
HEAR THE WORLD

 Sonova AG · Laubisrütistrasse 28
CH-8712 Stäfa · Switzerland

